



**Immune Technology Corp.**

*The Resource for Virology Research*

**SIVmac239 gp120  
ELISA Development Kit**  
Catalog Number: IT-E3Ag-SIVmac239gp120

**Description:** SIVmac239 gp120 ELISA Development Kit contains the key components required for the quantitative analysis of gp120(SIVmac239) concentrations in cell culture supernatants and serum within the range of 4-256 ng/ml in a sandwich ELISA format. The components supplied in this kit are sufficient to assay gp120(SIVmac239) in five 96-well ELISA plates.

**REAGENTS PROVIDED**

**Capture Antibody:** 100 $\mu$ l of 1mg/ml anti-gp120(SIVmac239) monoclonal antibody.

**SIVmac239 gp120 Standard:** 50 $\mu$ l of 50 $\mu$ g/ml recombinant gp120(SIVmac239).

**Detection Antibody:** 50 $\mu$ l of anti-gp120(SIVmac239) polyclonal antibody.

**Anti-Rabbit-IgG-HRP Conjugate:** 50 $\mu$ l of HRP-conjugated anti-rabbit IgG.

**RECOMMENDED MATERIALS & SOLUTIONS\***

ELISA 96-well plates (Corning Prod # 3590 or equivalents)

**Block Buffer:** 5% milk in PBS

**Wash Buffer:** 0.05% Tween-20 in PBS

**Diluent:** 0.05% Tween-20, 0.5% milk in PBS

**Substrate:** TMB Peroxidase Substrate

**Stop Solution:** 2N Sulfuric Acid


*\*Alternatively, these could be purchased under Cat.# IT-200-002 — ELISA Plate/Buffer/Substrate Kit.*

**PLATE PREPARATION**

1. For each 96-well plate, dilute 20 $\mu$ l of Capture Antibody with 10.5ml of 1xPBS to prepare a coating solution. Immediately add 100 $\mu$ l of the coating solution to each well. Seal the plate and incubate overnight at 4°C.
2. Remove the coating solution by aspirating or decanting. Invert the plate and blot it briefly against clean paper towels.
3. Add 300 $\mu$ l of Block Buffer to each well. Incubate for at least 1 hour at room temperature.
4. Aspirate to remove Block Buffer and wash the plate 4 times with 300 $\mu$ l of wash buffer per well.

**ASSAY PROCEDURE**

1. **Standard/Sample:** Dilute standard with Diluent to eight concentrations (256ng/ml, 128ng/ml, 64ng/ml, 32ng/ml, 16ng/ml, 8ng/ml, 4ng/ml, and 0ng/ml). Immediately add 100 $\mu$ l of Standard and sample to each well in triplicate. Incubate at room temperature for at least 1 hour.
2. **Detection:** Aspirate and wash plate 4 times. Dilute 10 $\mu$ l of Detection Antibody with 10.5ml of Diluent to prepare a detection solution. Add 100 $\mu$ l of the detection solution into each well. Incubate at room temperature for at least 1 hour.
3. **Anti-Rabbit-IgG Peroxidase:** Aspirate and wash plate 4 times. Dilute 10 $\mu$ l of anti-rabbit-IgG-HRP conjugate with 10.5ml of Diluent. Add 100 $\mu$ l into each well. Incubate at room temperature for 1 hour.
4. **Substrate/Stop:** Aspirate and wash plate 4 times. Add 100 $\mu$ l of TMB Peroxidase Substrate into each well. Incubate at room temperature for 20 minutes. Then add 100 $\mu$ l of Stop Solution to each well.
5. **Read:** Determine the optical density of each well within 30 minutes using a microplate reader set to 450nm.
6. **Analysis:** Average the triplicate reading for each standard, control, and sample, then subtract the average zero standard optical density. Create a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) or other curve-fit. The SIVmac239 gp120 concentration in sample can be determined by regression analysis. If samples have been diluted, the concentration calculated from the standard curve must be multiplied by the dilution factor.

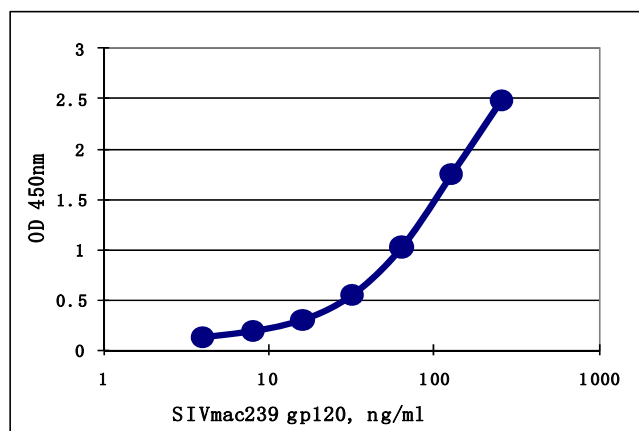

**Reference**

John R. Crowther. The ELISA Guidebook (Methods in Molecular Biology), Humana Press, 2000.